

ULTRA LOW TEMPERATURE FREEZERS -90°C Such Cryo-Pride





ilShinBioBase





0 🔤





DEEP FREEZER JABA SERIES

[JABA] Series

New JABA series of Cryo-Pride is equipped with cutting edge COSCOM (Combination Of Single Compressor Operation Mode) technology by two independent compressor mechanisms.

Combination of VIP(Vacuum Insulation Panel) and insulation by high pressure injection enable maximum insulation as well as lowest energy consumption for keeping your valuable samples with maximum safety. Auto Cascade system with COSCOM technology protects your samples from any risk other than blackout. With low noise and all time monitoring "Blackbox" service, our JABA series shall be successful companion of your comfortable research.




Main Features

JABA SERIES

No worries, No hassles in any emergency-**1. INDEPENDENT COOLING SYSTEM**

COSCOM System

Two independent cooling mechanisms. One cooling systems maintains -80°C in case of system failure.

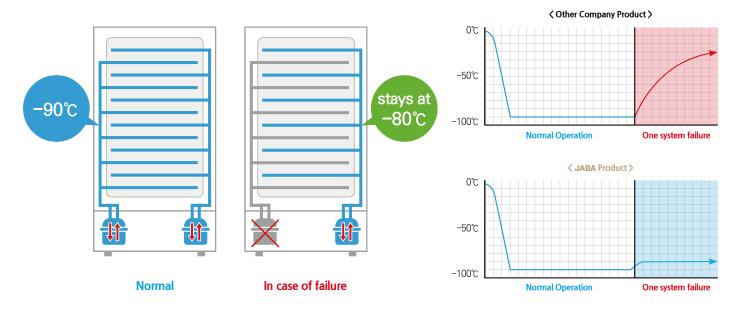

2. 24/7 real time monitoring **BLACK BOX** service

GLOBAL REAL TIME CARE SERVICE






Remote monitoring by Internet of things technology

Manufacturer's central monitoring server gives 24/7 real time monitoring by IoT technology. In case of any event, server gives notices to user and designated service team.

Central server monitoring data

- 1. Temperature: Chamber temperature of the freezer
- 2. Trend information: Temperature variation for last 24 hours
- 3. Alarm event: Compressor over temperature, compressor over current, system overheat, temperature sensor failure, power off, low voltage, cold alarm, warm alarm and door open alarm.

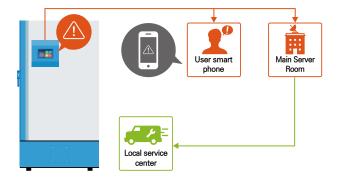

Prompt and proactive action

User can find any events by application of smart phone. Diagnosis and corrective action are noticed.

*App. OS: Android 2.3 (Gingerbread) or above, IOS 9.2 or above

Any event alarm goes to user' smart phone and to main server. Main server shall review history data of the unit and prompt diagnosis is executed. Service action is instructed proactively!

< Conventional service protocol >

Step	Processing details
1	Checking symptoms
2	Service request
3	Engineer field visit and diagnosis
4	Service parts request
5	Engineer visit for repair
6	System reactivated

Complication, time consumption

< New JABA service protocol >

Step	Processing details			
1	Event happens, service instruction provided			
6	Completion of service			

No hassle, Non-Stop service!

No need for searching your samples!

1. INVENTORY MANAGING SYSTEM

Inventory system

Digital memo function for sample storage location.Energy saving with minimized door open time.



4. data review by touch and scroll- TREND DATA



🔹 Graph

- Review on the screen over 43,000 hours past trend.
- USB download of data is available.





5. monitoring in your hand- MOBILE CONTROL

Mobile remote control

Monitoring and control by mobile application is available with your computer, Lap top, smart phone or tablet PC.




CONTROL SYSTEM MOBILE APPLICATION



* Main

* Log in Mobile App account open and


log in.

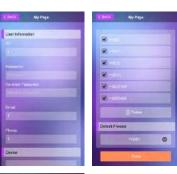

* My page

Registration of user unit. Service authorization.

Setup, My page, Temperature graph, Event menu. Real time monitoring page.

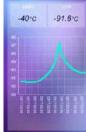

* History

Present data and trend data for last 43,000 hours on screen.

* Event Alarm events display and review.



* Setting

Parameter setup, password setup and authorization control.

Better protection and lower energy consumption! **POWERFUL INSULATION**



Improved insulation

Combination of VIP and PU provides better insulation that enables fast recovery of chamber temperature.

-50°C

-100°C

Normal Operation



slow recovery, high energy consumption, risk on samples

< Other Company


Product >


Temperature change before and after door open/close

Temperature variation during power off






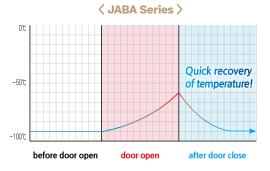

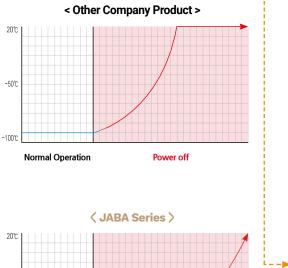

(about 11 hrs)

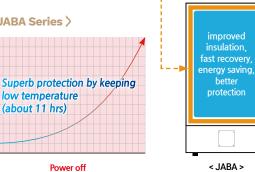

Don't get bothered by loud noise anymore **LOW NOISE STRUCTURE**



Designed for low noise

Comfortable work environment is provided. Noise level is as low as home refrigerator.

Low noise below 55.1dB(A)



< TUV NORD Performance Certificate >

DEEP FREEZER GUDERO PLUS SERIES

[GUDERO PLUS] Series

Gudero Series was introduced in 1988 and has been upgraded since then. New Gudero Plus Series is equipped with COSCOM technology for double safety and robust quality. Easy user interface is beneficial for long term depository by simplified control parameters with only 3 buttons.



Main Features

GUDERO PLUS SERIES



Independent cooling mechanisms

1. DOUBLE SAFETY

*COSCOM system

Even with failure of one compressor, temperature is still maintained at -80 $^\circ\!C.$

Easy Control by 3 Buttons

2. SIMPLIFIED CONTROL SYSTEM

*Micro process control system

Gudero Plus is optimal for long term bio-depository with easy user interface and simplified performance.

No worries about overheat or dust accumulation! 3. PIN-LESS CONDENSATION SYSTEM

*Pin-less

No need for cleaning, no need for filter replacement. Compressor overheat is prevented.

DEEP FREEZER GUDERO CHEST SERIES

[GUDERO CHEST] Series

Gudero Chest Series is optimal for long term storage of large size samples at the range of -90 $^{\circ}$ -120 $^{\circ}$.

-10 °C \sim -120 °C Gudero Chest lineup is with reliable cascade refrigeration technology. With accumulated know-how, micro process control system provides easy and hassle free long term bio-depository solution.



Main Features

GUDERO CHEST SERIES

Cryo-Pride

Light and handy

. POLYETHYLENE INNER DOOR

Light and handy inner doors make your research more comfortable.

Self-diagnosis function, smart & easy maintenance

2. SELF-DIAGNOSIS FUNCTION

Temperature and voltage parameters are monitored and alarm is generated at any emergency situation.

Temperature check by one glance MONITORING SYSTEM

RS-232 and 485 communication port allows real time monitoring by personal computer.

GUDERO - CHEST





Main Features

POCKET SERIES



Convenient access

. COMPACT DESIGN

Table top or under counter storage is possible even with small space. Fast and easy access is privilege of Pocket series.

If you need more storage?

2. 2nd STACKING

Easy build up without need of more space.

QUICK FREEZER

[QUICK FREEZER]

General use of deep freezer is long term storage of samples but special samples such as blood and specific reagents need prompt freezing for keeping their properties.

QDF Series is optimal solution for fast and reliable freezing without chance of any damage on your samples.



Main Features

QUICK FREEZER

Powerful & fast freezing system **1. FAST FREEZING SYSTEM**

Specially designed refrigeration system and powerful ventilation pull down from 20°C to -40°C within 1 hour.

Perfect insulation 3. 4 LAYER GASKET

00

3 layers on chamber side and 1 on door side, minimizing loss of cold air.

Rail type air hole drawer with enough space **2. REMOVABLE RAIL DRAWER**

Rail type drawer with air hole enables good circulation of cold air. Convenient way of categorized sample storage.

DETAIL & FEATURES DEEP FREEZER



Smart Monitoring Application

24/7 monitoring available on your smart phone. Proactive notice from command room in case of alarm event.



Finger Print Door Lock

Enhanced security: by finger print registration or by card key.



Inner door for each shelf

HDPE inner doors minimize cold air loss.



Powerful insulation

High pressure injected POLY URETHANE with 130mm thickness keeps temperature less than 3°C between chamber inside and outer wall.





Control panel

 TFT LCD 7" Touch Screen
 Smart Microprocessor control



Sliding Rack

Convenient storage by separated racks with sliding mode.



Anti-explosion door lock

Auto relief system(open & close) works in case of high pressure formed inside of chamber. Key set is provided for locking function.


Black Box

Manufacturer's central monitoring server gives 24/7 real time monitoring by IoT technology. In case of any event, server gives notices to user and designated service team.

< CENTRAL SERVER MONITORING DATA >

1. Temperature: chamber temperature of the freezer.

Trend information: temperature variation for last 24 hours.
 Alarm event: compressor over temperature, compressor over current, system overheat, temperature sensor failure, power off, low voltage, cold alarm, warm alarm and door open alarm.

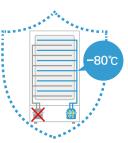

Coscom system (Combination of single compressor operation mode)

Even with failure of one compressor, temperature is still maintained at -80 °C. Fast pull down to -90°C and energy efficient ECO mode.

SPECIFICATIONS DEEP FREEZER SERIES

JABA Series

Series	DF8510	DF8514	DF8517	DF8520	DF8524	DF8530	
Capacity	276 liters	369 liters	484 liters	605 liters	730 liters	832 liters	
Ext. Dimensions (WxDxH mm)	920 x 805 x 1660	920 x 805 x 1980	920 x 955 x 1980	1070 x 955 x 1980	1227 x 957 x 1980	1227 x 985 x 1980	
Int. Dimensions (WxDxH mm)	600 × 480 × 970	600 x 480 x 1290	600 x 630 x 1290	750 x 630 x 1290	890 x 635 x 1290	890 x 725 x 1290	
Temperature Range (℃)	⊃°09- ~ ⊃°03-						
Insulation		130 m	m Non CFC, foamed ir	n place polyurethane w	rith VIP		
Electrical			110V, 220V 60	Hz, 230V 50Hz			
Certification		ISOS	9001 / ISO14001 / ISO)13485 / CE / KGMP /	ETL		
Weight (kg)	256	279	307	340	378	410	
Refrigeration	2 x 1.25 HP						
2" Cryo Boxes	15 boxes x 12 racks	15 boxes x 16 racks	20 boxes x 16 racks	20 boxes x 20 racks	20 boxes x 24 racks	25 boxes x 24 racks	
3" Cryo Boxes	9 boxes x 12 racks	9 boxes x 16 racks	12 boxes x 16 racks	12 boxes x 20 racks	12 boxes x 24 racks	15 boxes x 24 racks	
2" Boxes/Freezer	180 boxes	240 boxes	320 boxes	400 boxes	480 boxes	600 boxes	

*Voltage range: 110V(100-120), 220V(208-230), 230V(220-240)

GUDERO PLUS-UPRIGHT Series

Series	DF3510	DF3514	DF3517	DF3520	DF3524	DF3530
Capacity	276 liters	369 liters	484 liters	605 liters	730 liters	832 liters
Ext. Dimensions (WxDxH mm)	920 x 805 x 1660	920 x 805 x 1980	920 x 955 x 1980	1070 x 955 x 1980	1227 x 957x 1980	1227 x 985 x 1980
Int. Dimensions (WxDxH mm)	600 × 480 × 970	600 × 480 × 1290	600 × 630 × 1290	750 x 630 x 1290	890 x 635 x 1290	890 x 725 x 1290
Temperature Range (℃)	-10°C ~ -35°C					
Insulation		130 m	m Non CFC, foamed ir	n place polyurethane w	vith VIP	
Electrical			110V, 220V 60	Hz, 230V 50Hz		
Certification		19	SO9001 / ISO14001 /	ISO13485 / CE / KGM	IP	
Weight (kg)	212	233	260	291	326	360
Refrigeration	1 x 1.25HP					
2" Cryo Boxes	15 boxes x 12 racks	15 boxes x 16 racks	20 boxes x 16 racks	20 boxes x 20 racks	20 boxes x 24 racks	25 boxes x 24 racks
3" Cryo Boxes	9 boxes x 12 racks	9 boxes x 16 racks	12 boxes x 16 racks	12 boxes x 20 racks	12 boxes x 24 racks	15 boxes x 24 racks
2" Boxes/Freezer	180 boxes	240 boxes	320 boxes	400 boxes	480 boxes	600 boxes

GUDERO PLUS-UPRIGHT Series

Series	DF8510	DF8514	DF8517	DF8520	DF8524	DF8530
Capacity	276 liters	369 liters	484 liters	605 liters	730 liters	832 liters
Ext. Dimensions (WxDxH mm)	920 x 805 x 1660	920 x 805 x 1980	920 x 955 x 1980	1070 x 955 x 1980	1227 x 957 x 1980	1227 x 985 x 1980
Int. Dimensions (WxDxH mm)	600 × 480 × 970	600 × 480 × 1290	600 x 630 x 1290	750 x 630 x 1290	890 x 635 x 1290	890 x 725 x 1290
Temperature Range (℃)	-2°06- ~ 3°09-					
Insulation		130 m	m Non CFC, foamed ir	n place polyurethane w	vith VIP	
Electrical			110V, 220V 60	Hz, 230V 50Hz		
Certification		ISOS	9001 / ISO14001 / ISO)13485 / CE / KGMP /	' ETL	
Weight (kg)	240	262	289	320	355	390
Refrigeration	2 x 1.25HP (1 x 1.25HP)					
2" Cryo Boxes	15 boxes x 12 racks	15 boxes x 16 racks	20 boxes x 16 racks	20 boxes x 20 racks	20 boxes x 24 racks	25 boxes x 24 racks
3" Cryo Boxes	9 boxes x 12 racks	9 boxes x 16 racks	12 boxes x 16 racks	12 boxes x 20 racks	12 boxes x 24 racks	15 boxes x 24 racks
2" Boxes/Freezer	180 boxes	240 boxes	320 boxes	400 boxes	480 boxes	600 boxes

GUDERO-CHEST Series

Series	DF4510	DF4517	DF9010	DF9017	DF12010
Capacity	295 liters	295 liters	295 liters	492 liters	295 liters
Ext. Dimensions (WxDxH mm)	1590 x 790 x 1030	2180 x 790 x 1030	1590 x 790 x 1030	2180 x 790 x 1030	1590 x 790 x 1030
Int. Dimensions (WxDxH mm)	870 × 500 × 670	1470 × 500 × 670	870 × 500 × 670	1470 × 500 × 670	870 × 500 × 670
Temperature Range (℃)	−10°C ~ −45°C		-40℃ ~ -90℃		−80℃ ~ −120℃
Insulation		130 mm Non C	FC, foamed in place polyur	ethane with VIP	
Electrical		1	10V , 220V 60Hz , 230V 50	Hz	
Certification	ISO9001 / ISO 140	01 / ISO 13485 / CE	ISO9001 / ISO14001 / ISO13485 / CE / KGMP		
Weight (kg)	207	257	236	286	280
Refrigeration	1 x 1.25HP		2 x 1.25HP		2 x 1.5HP
2" Cryo Boxes	12 boxes x 18 racks	12 boxes x 30 racks	12 boxes x 18 racks	12 boxes x 30 racks	12 boxes x 18 racks
3" Cryo Boxes	8 boxes x 18 racks	8 boxes x 30 racks	8 boxes x 18 racks	8 boxes x 30 racks	8 boxes x 18 racks
2" Boxes/Freezer	216 boxes	360 boxes	216 boxes	360 boxes	216 boxes

POCKET Series


Series	DF8502	DF8503		
Capacity	40 liters	80 liters		
Ext. Dimensions (WxDxH mm)	470 × 620 × 780	770 × 620 × 780		
Int. Dimensions (WxDxH mm)	300 × 440 × 300	600 × 440 × 300		
Temperature Range (℃)	−60°C ~ −86°C			
Insulation	High Density poly	urethane with VIP		
Electrical	110V, 220V 60	110V, 220V 60Hz, 230V 50Hz		
Certification	ISO9001 / 14001 /	/ 13485 / CE / ETL		
Weight (kg)	79	99		
Refrigeration	3/8 HP	5/8 HP		
2" Cryo Boxes	15 box x 2 racks	15 box x 4 racks		
3" Cryo Boxes	9 box x 2 racks	9 box x 4 racks		
2" Boxes/Freezer	30 boxes	60 boxes		

QUICK FREEZER

VOICH FREEZER				
Series	QDF4517			
Capacity	484 liters			
Ext. Dimensions (WxDxH mm)	920 x 955 x 1980			
Int. Dimensions (WxDxH mm)	600 × 630 × 1290			
Temperature Range (℃)	−10°C ~ −40°C			
Insulation	130 mm Non CFC, foamed in place polyurethane with VIP			
Electrical	380V, 3Ph, 60Hz			
Certification	ISO9001 / 14001 / 13485 / KGMP			
Weight (kg)	307			
Refrigeration	5 HP			
Drawers/Freezer	total 11 Drawers			

*Voltage range: 220/380V, 3Ph, 50/60Hz

DIMENSIONS & STORAGE CAPACITY Preview

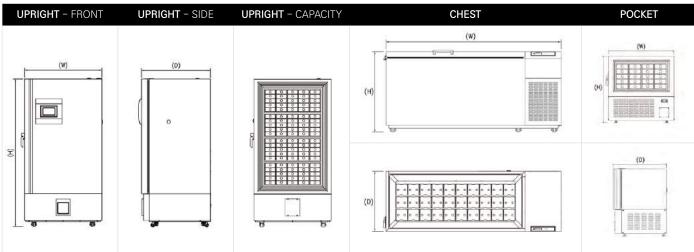

*Stainless steel option: All models have the option for a stainless steel interior. Please add "S" behind the model number to order stainless steel interior, (for example DF 8524S)

CONFIGURATION

for FREEZER SYSTEM



environment

ACCESSORIES

for FREEZER SYSTEM

Rack Series


4

8

11

Item NO.	Boxes Size	Quantity	Purpose
CR0060	2"	12	for Chest Freezer
CR0061	3"	8	TOI Chest Theezer
SR0062	2"	15	Upright Freezer
SR0063	3"	9	xx02, xx03, xx10, xx14
SR0064	2"	20	Upright Freezer
SR0065	3"	12	xx17, xx20. xx24
SR0066	2"	25	Upright Freezer w/20
SR0067	3"	15	Upright Freezer xx30

Purpose	Quantity	Boxes Size	Item NO.	
for Chest Freezer	12	2"	CR0060	
TOI CHEST HEEZEI	8	3"	CR0061	
Upright Freezer	15	2"	SR0062	
xx02, xx03, xx10, xx14	9	3"	SR0063	
Upright Freezer	20	2"	SR0064	
xx17, xx20. xx24	12	3"	SR0065	
Upright Freezer xx30	25	2"	SR0066	
Oprigrit i reezer xx30	15	3"	SR0067	

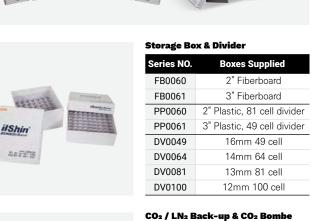



02 / LN2 Back-up & CO2 Bombe					
Series NO.	Description				
LC0070	CO ₂ back-up system with battery and temperature display				
LN0070	LN ₂ back-up system for protection of samples from emergency situation (black out, temperature drop)				
CO2000	Carbonic Acid gas tank 20KG/CAN: Siphon included, Unit: 20KG/CAN				

Description

Internet based real time

Proactive service for any

monitoring (24/7).

emergency.

Black Box Terminal Unit

Series NO.

BBD1000

10



Series NO. Description

Control System

Security Finger Print & Card Access

Enhanced security by finger-SCA0002 print and access card key to improve D.A.S function.

Safety Door Lock Device

Dialing System

Series NO. Description Auto relief system(open & close) works in case of high pressure (above 0.0008kg/ SD4989 cm²) formed inside of chamber. Key set is provided for locking function.

Series NO. Description Wireless sensor based emergency protection by DS0119 giving alarms up to 5 points telephone numbers or messages.

	AVR Device	•		
escription	Series NO.	Description		
voltage under ge environment	AVR0003	Regulating voltage supply under unstable voltage		

12

AVR0005

Chart Reco	Safety G	
Series NO.	Description	Series N
CH0067	Temperature recorder. 24/7 day circular chart recorder	GL0005

Safety Glove				
Series NO.	NO. Description			
GL0005	Protection from frost bite			

Voltage Booster						
Series NO.	Description					
VB0081	Boosting voltage under low voltage environment (AC15V)					





动车扩新部

3

7

9





Product Inquiry

Base	Address	Telephone	Fax
Overseas Sales Department	Lee Jong Ho / Deputy General Manager – jhlee@1sbb.com 84 Samyuksa-ro 548 Street, Dongduchon city, Kyungki do, Korea. ZIP 11341	+82-70-4354-3951	+82-70-7950-3951

ilShinBioBase Locations

Base	Address	Telephone	Fax
본사 및 생산	경기도 동두천시 삼육사로 548번길 84	031-867-1384	070-7950-3911
서울/경기지사	서울특별시 중랑구 중랑천로 77	1577-4053	02-491-4073
수원지사	경기도 화성시 봉담읍 안당하리길 34	031-298-8147	031-298-8149
대전지사	대전시 유성구 테크노8로 58 1F	042-824-1145/6	042-824-1147
영남지사	대구광역시 동구 동호로 75, 4F	070-4354-3977/5000	070-7950-3941
호남지사	광주광역시 광산구 신가삼효로 20-15	062-951-8010	062-951-8011